


SSA-806 Installation

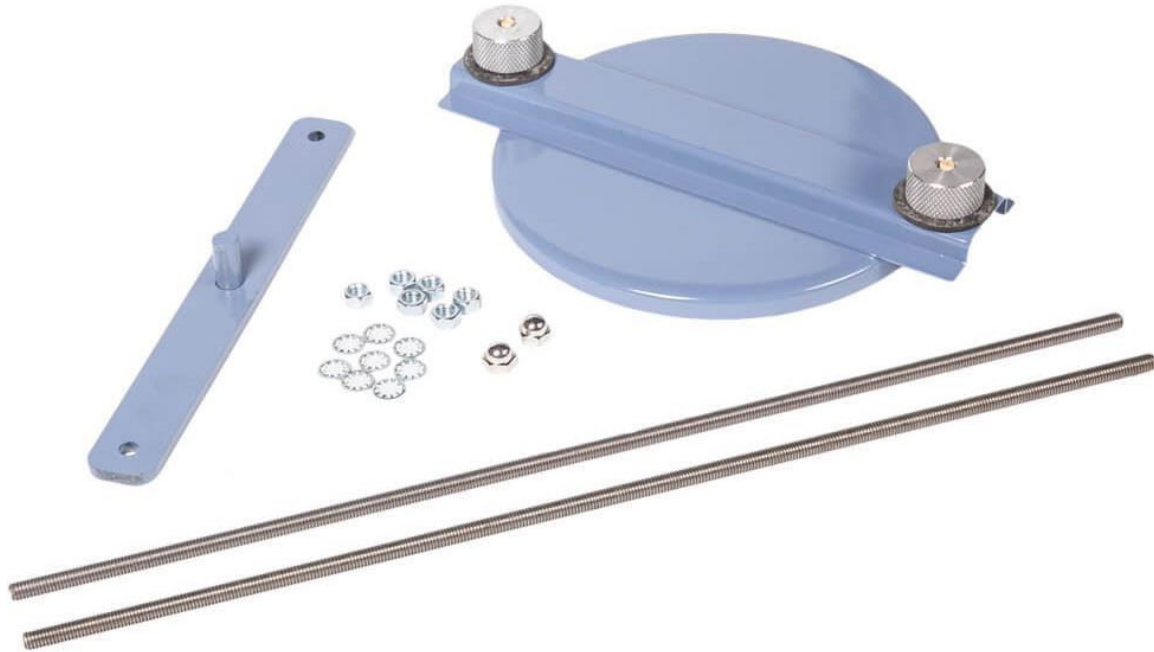

SSA-806

INSTALLATION

Remove existing clamp rods and cover from your existing SS-15 or SS-15D.

Replace the clamp rods with full threaded rods, place lock washers and hex nuts on rods and turn them on the rods about $\frac{1}{2}$ " from the ends of the threaded rods. Place threaded rods on the platform, thread lock washer and hex nut on rods below platform and tighten all nuts to hold clamp rods securely in place.

Take EZ Clamp top cover and push the buttons on the EZ Clamp knobs, allowing the spacers to drop out of the EZ Clamp knobs. Still holding the buttons, place threaded clamp rods through the EZ Clamp knobs and slide to desired position and release buttons. The cover should remain in place. Place lock washer

and hex nut on top of rods about $\frac{1}{2}$ " from end. Place top support bar on rods. Place acorn nuts on top of threaded clamp rods and tighten all nuts.

To move EZ Clamp cover, just push the buttons on EZ Clamp knobs and slide up and down.

Place stack of sieves on the platform and simply push the EZ Clamp knob buttons and slide cover down on top of the sieve stack and release. Finally, turn each EZ Clamp knob until the EZ Clamp knob tightens to hold the sieve stack.

Release by unscrewing the knobs enough to release tension on the stack, the push the buttons on the EZ Clamp knobs and slide the cover up to desired position and remove sieve stack.

